

Step 1 Check required specification for Phone or Tablet

① Before installing OKOK App, kindly check if the phone or tablet meets following requirement:



Compatible with iOS7.0 and above with Bluetooth 4.0 and above

Phone 4s, iPhone 5, iPhone 5c, iPhone 5s
iPhone 6, iPhone 6 Plus, iPod touch 5, iPad 3
iPad 4, iPad Air, iPad Air 2, iPad mini
iPad mini 2, iPad mini 3




Compatible with Android 4.3 and above with Bluetooth 4.0 and above

SAMSUNG: galaxy s5, galaxy s4, galaxy s3, galaxy note4, galaxy note3, galaxy note2, galaxy grand2
XIAOMI: XIAOMI Note, XIAOMI4, XIAOMI3, XIAOMI2/2s, XIAOMI tablet, HONGMI Note 4G, HONGMI 1s
HUAWEI: Ascend P8, Ascend P7, Honor 6Plus, Honor 6, mate7, Ascend G7

Support Android4.3 and above for selected Android devices with Bluetooth 4.0 and above

Step 2 Install the batteries

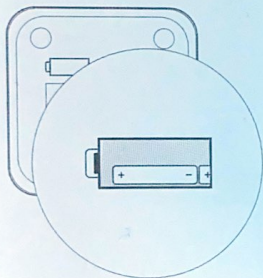

2

- 1/ - Open the cover of battery
- 2/ - Install the batteries in the right direction
- 3/ - Close the cover

Step 3 Install App

- 1/ Switch on Bluetooth


Use your device to scan QR Code or search OKOK International from App Store/Google Play to install App.

- 2/ Phone or Tablet connects to Net (WiFi, 3G or 4G)




Step 4 Pair the devices, register and complete user information


Open APP, begin to weigh and pair

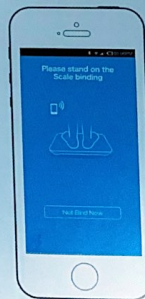

Register as a user

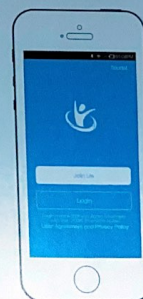

Complete personal information

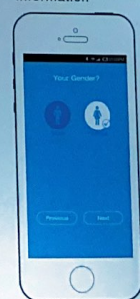


3

Step 5 Weighing

Enter into "Dynamic Weight" interface



Stand on the scale and synchronize weight data onto APP

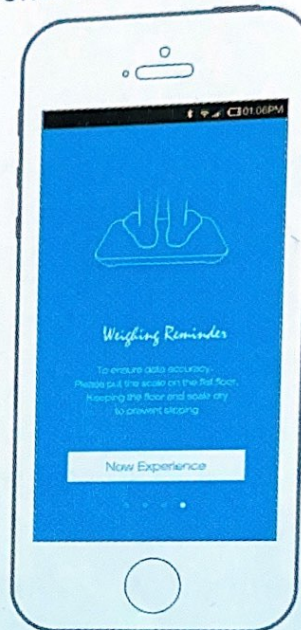

Step 6 View the test result

Enter into "Weight Trend" interface



View test report

