


Product function

1. USB contact charging:
Charging time is less than 120 minutes, once charged, it can last for more than one week, saving energy and environmental friendly.
2. convenient to use:
Simple operation, four mode options: A, strong brushing B, comfortable brushing C, high-frequency massage D, whitening teeth
3. brush your teeth in all directions:
After the switch is activated, it takes only 45 seconds to complete a 360-degree brushing action. (45 seconds automatic stop)
4. automatic memory function:
Intelligent store user habits - efficient and convenient.
5. Three in one functions of brushing, massage and whitening teeth.

Product Features

- 1 this product adopts all-round waterproof seal design, LPX7-convenient, durable and clean
- 2 U-shaped toothbrush head design, using food-grade silicone antibacterial material, 360 degrees contacting with gums and teeth. Brushing your teeth and brings a quick and comfortable experience to your mouth every time.
- 3 the use of this product can prevent gingival recession, promote periodontal blood circulation and reduce the occurrence of periodontal disease.

Introduction of product accessories

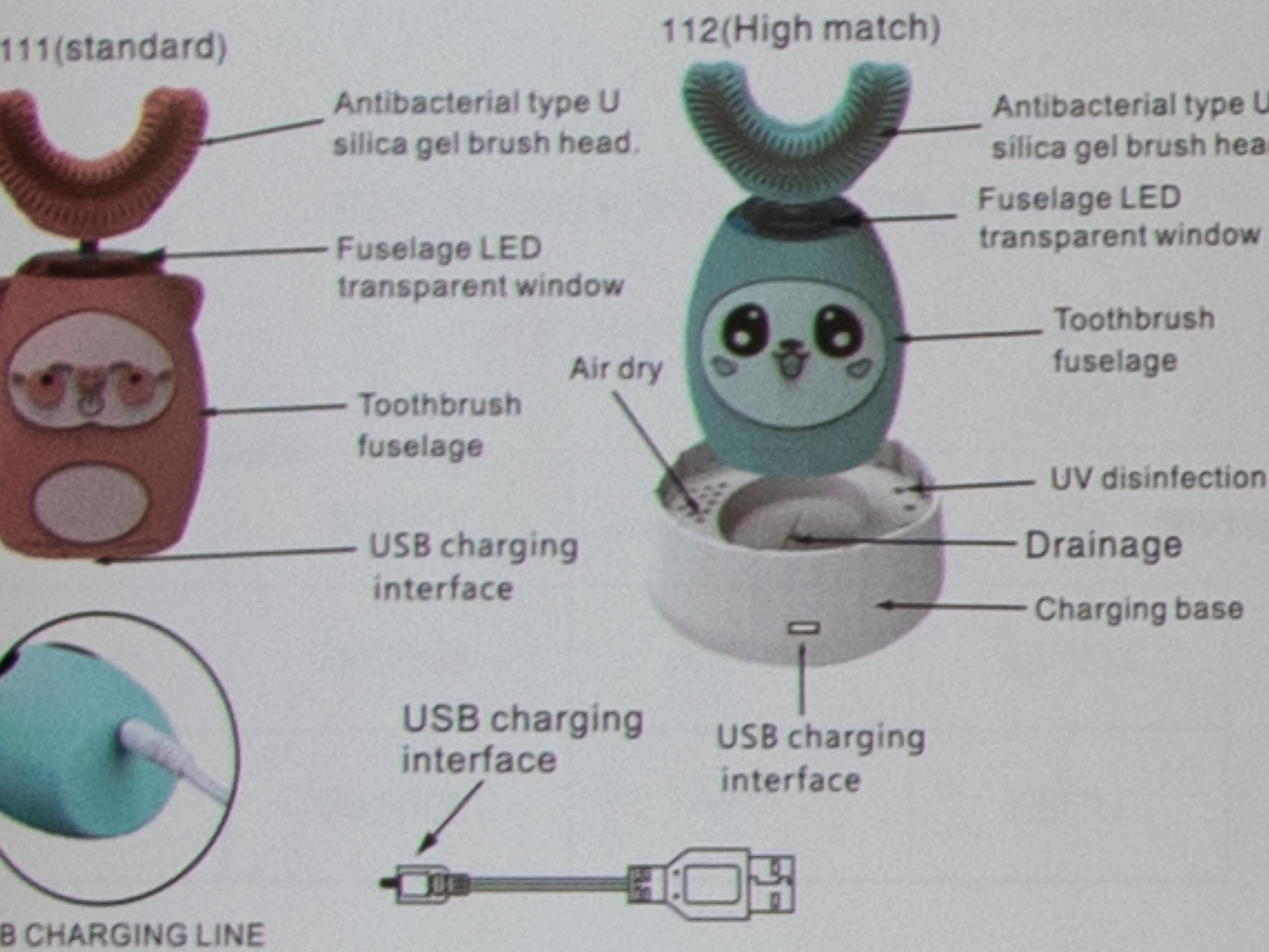

2. Charging description

First, a small USB at one end of the USB data line is inserted into the plug socket of the charging base, and the other end is connected to the USB power supply. Then the mainframe is placed vertically on the charging base grooves to recharge.

- ① when the fuselage electrode is in good contact with the charging base, the charging base indicator light is lit. When the fuselage led is not bright or blink it is indicated that the fuselage is not placed well, please replace it.
- ② Charging, the fuselage LED glittering, full of light.

3. Product parameters

Product size	108x45x33mm	Product weight	92g±10g
Rated voltage	3.2~3.7V	Rated power	<5W
Working current	<1.6A	Charge time	<120Minute

3. Use steps

1. Parallel to the body and align the hole with the power shaft, and put it in place.

2. This product can use liquid toothpaste or ordinary toothpaste.

3. When using ordinary toothpaste, you need to wet the braces, then put them into the mouth, let the upper and lower teeth gently fit on the alveolar, and turn on the switch, shake the hands to touch the left and right teeth, the best frequency is 63 times / minute, A comfortable and fast tooth cleaning process can be completed in 45 seconds. The braces can be modified to suit the length.

4. Clean the braces and the body, braces with water rinse, use dry cloth to wipe the fuselage.



The mouthpiece can be modified to a suitable length according to your mouth.

Matters needing attention

- 1 the normal working temperature of the product is 0 ° C -50 ° C.
- 2 when the product is charging, due to the characteristics of the product, please charge at room temperature.
- 3 if the product is not used for a long time, it is recommended to charge and discharge once within three months.
- 4 the product should be stored in a dry, cool place, avoiding the sun.
5. After brushing your teeth, please use a damp cloth to dry the water and stagnant water on the body, and try to avoid be emerged in the water.
- 6 the product has large output power, do not work for a long time. Usually less than 10 minutes.
7. Please keep the body dry when charging. The power supply is 5v0.5a to 5v1a charging power supply.
- 8 the picture is for reference only, please base on physical product

Product warranty card

NO.:	
Date of purchase:	
Customer:	
Dealer:	
Name:	Address:
Telephone:	Name and address:
Batch Number:	

Product Name: Intelligent U-shaped Electric Toothbrush For Children
Rated voltage: 5V
Rated power: 5W