


Smart Neck Massager


HX-1680

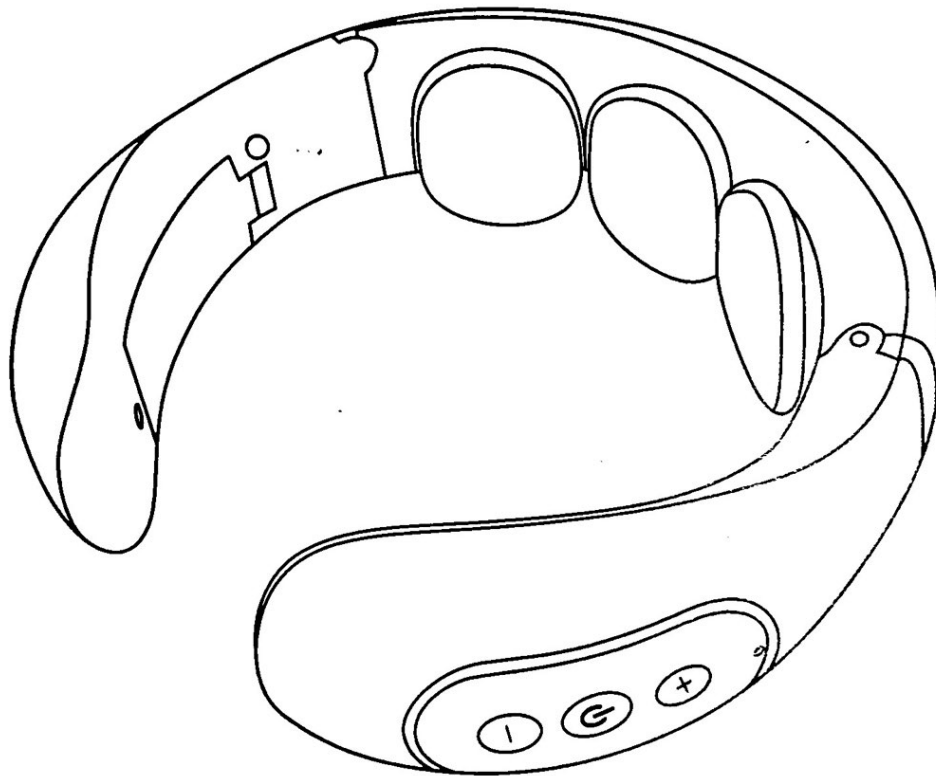

an instruction manual

Product features and functional features

1. According to ergonomics design, beautiful and scientific, in line with the neck curve;
2. Built in rechargeable lithium battery for power supply, safe, durable and reliable;
3. Adopt 3D intelligent fitting technology, automatically adjust massage head according to neck curve, or choose to use External electrode piece massages other parts;
4. 4 kinds of massage techniques, 15 levels of intensity, can be selected according to different needs, suitable for a wider range of people;
5. Massage the acupoints with external electrodes;
6. The comfort of hot compress function was further improved.

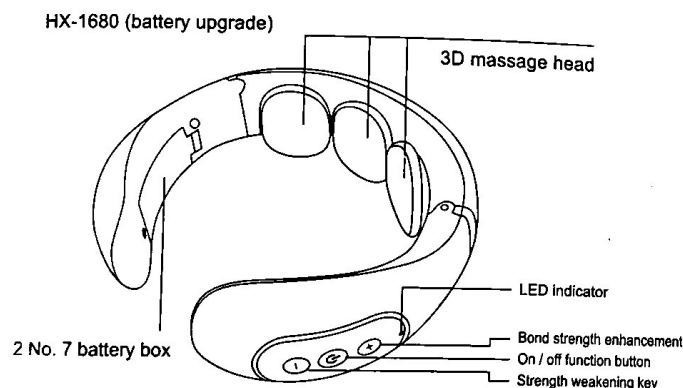

Operation method of hx-1680 (battery upgrade) function keys

How to turn it on

Open the cover of the battery box, and put in two No.7 batteries according to the marked positive and negative poles in the battery box. Pay attention to the electricityCell polarity to prevent incorrect installation, and then close the battery box cover.

※ Long press the @ key on the instrument for 2 seconds, the instrument starts to work; the LED indicator light is on, and the default kneading mode is on.

※ Mode selection: short press the @ key to select the corresponding working mode.

The modes are hammering mode and acupuncture mode. Mode, automatic mode, kneading mode, (the automatic mode is: the above three modes of automatic cycle.)

Strength selection

※ ⊕ The bond is slightly increased, which can increase the strength;

※ ⊖ The bond is the strength weakening key, which has the weakening strength about.

※ In the process of use, after selecting the working mode, press the ⊕ key once to increase the intensity of one gear, a total of 15 gears

Strength: press the ⊖ key once to reduce the intensity of a gear, please choose the comfortable strength you accept.

How to shut down

※ Press and hold the @ key on the instrument for 2 seconds to turn off the machine and the machine will stop working. (if it is not used for a long time, please take out the No.7 battery)

Open the cover of the battery box, take out the two No. 7 batteries, and then close the cover of the battery box.

Product model and parameters

Battery upgrade
Product model: hx-1680
Input power supply: 3V=12 sections, No.7, 1.5V
Alkaline battery
Rated power: 4W
Timing: 15 minutes
Net weight: 140g

Preparation before use

1. Clean the neck hair and clothes, otherwise it will affect the comfort of the product;
2. When massaging other parts of the body with electrodes, please stick the two electrodes directly to the skin at the same time. When you need to adjust the massage position and electrode piece, please turn off the instrument and turn it on again after adjustment. Start operation.

How to operate



reminder:

Because each person's skin nature is different, the effect of using will also be different, dry skin. The resistance is small, if the force is too weak or tingling during normal use, it is recommended to use warm towel to compress the massage part for about 1 minute! (oily. This is not necessary for neutral skin)

Product introduction

The wireless intelligent cervical massage device adopts the streamline circular design which conforms to the physiological curvature of human cervical vertebra. Rub neck and shoulder. To deal with neck acid, numbness, distension, stiffness. Promote neck, shoulder and other parts of the muscle relaxation. During the whole process, keep the correct physiological posture all the time. Fit, massage more comfortable. Wider neck circumference adjustment angle. It can be used by the whole family. It is fit for people. Increase shoulder massage output, massage part expansion, electrode can massage arm, waist, leg and other parts. Light and small. Smart, easy to carry, massage experience anytime and anywhere, travel and home are available. Multiple massage modes can be switched with one key. For neck relaxation massage care, one machine multi-purpose. According to the neck physiological curve and elastic bending design, press the With high comfort and wide applicability, massage every day and take care of neck every moment.

Scope of application

Deep care of cervical vertebra, relax local muscles, to eliminate fatigue.

reminder

1. Before using the product, comb the neck, hair and jewelry, etc. to prevent the product from being blocked and unable to work normally;
2. After the product is turned on, the strength can be gradually adjusted by the "+" - " button according to the comfort degree.



Before using hx-1680 series products for the first time, it is better to charge the products for 2 hours.

matters needing attention

Usage environment

1. Do not use this product in the environment with excessive humidity or dust concentration, so as to avoid failure;
2. Please do not use this product in the indoor temperature above 40 °C, and it is not suitable to use it at minus 10 °C;
3. Please do not place the product near electric heater or stove or in the place where there is direct sunlight and water;
4. This product is designed to be used by people and can not be used for any commercial purpose;
5. If the product is moved from low temperature storage environment to warm environment, please put the product normal room first. Keep it in a warm environment for 1 hour, so as not to affect the normal use of the machine due to the water vapor condensed by the air. To cause machine failure;
6. The product has a heating surface, and people who are not sensitive to heat must pay attention to it. 7. Remove the battery from the appliance before abandoning the appliance; 8. When taking out the battery, the appliance must be powered off;
9. The battery should be placed safely. Throw it into the red non reusable garbage can. Do not throw it randomly to avoid environmental pollution.

DANGER

• Please never use it with the following medical electronic devices

1. Implantable medical electronic devices such as cardiac pacemaker;
 2. Life sustaining medical electronic devices such as artificial heart and lung
 3. Electrocardiograph and other medical electronic devices
- ※ It will cause malfunction of these medical electronic devices and may endanger life.
- If the patient uses the high-frequency surgical equipment and the product simultaneously, the electrode area of the product may lead to burns. Burns may also damage the product
 - If the product is used near shortwave or microwave therapy equipment, the output of the product may be unstable.
 - avoid contact with wounds and cracks during treatment
 - When using the product, the electrode should not be placed near the chest and heart, otherwise it will increase the risk of heart fibrosis.
 - The product is recommended for one user. Please clean the electrode on the product before use.
 - When using the product, the electrode sheet should be tightly attached to the skin and evenly contacted.
 - Stop using if you feel unwell during or after use and consult your doctor. Please consult your doctor for the following conditions. Otherwise,
 - Please consult your doctor for the following conditions. Otherwise, accidents or physical discomfort may occur.
1. Patients who are receiving medical treatment and feel abnormal.
 2. Patients with abnormal heart, brain and nerve.
 3. Patients with body temperature over 38 degrees (fever period).
 4. Patients with infectious diseases.
 5. Patients with skin perception disorder or abnormal skin.
 6. Patients who are taking drugs.
 7. Patients with abnormal blood pressure.
 8. Patients who are unable to express themselves freely.
 9. Patients with disturbance of consciousness caused by diabetes mellitus.
 10. Patients with skin allergy to stainless steel electrodes.
- It is forbidden to use in the following situations
1. Patients with malignant tumor
 2. Patients with bleeding tendency
 3. Pregnant women
 4. Patients with acute suppurative inflammation
 5. Pacemaker users
 6. Patients with metal implants
 7. Near the heart of the human body
- Do not use it near the heart, above the neck, head, mouth, pudendum or skin diseases, and cannot be used at the same time. Treat the soles of both feet. In addition, do not compress the viscera. Otherwise, accidents or physical discomfort may occur;
- Do not use with, other medical devices or more than 2 sets of this product, or it may cause nausea or discomfort.
- Do not use it for purposes other than massage recorded in the instruction manual, otherwise accidents may occur.

Precautions for using dry battery

1. Do not charge the battery, which may lead to the danger of battery leakage.
2. The rechargeable battery should be taken out of the appliance and put into the special charger for charging.
3. Pay attention to battery polarity to prevent wrong installation.
4. Do not short circuit or heat the battery. Throw into the fire and try to tear it apart.
5. Do not mix old and new batteries.

Use of electrode sheet and others

1. How to use electrode piece

※Connect the physiotherapy line with the socket of the physiotherapy line, and connect the other end with two pieces of electrode pieces, and stick the electrode piece close to the part to be massaged.

※Note: two electrodes must be pasted on one person before massage.

※Note: when the external electrode piece is used, the massage head of the instrument does not work.

※Attachment: massage electrode patch is consumable. Please purchase and replace it in time if it is used frequently for a long time.

2. In the process of massage, the neck has tingling pain and discharge inductance

※Due to different physical fitness, dry people have dry skin and low resistance.

Suggestion: wipe the neck skin with water or moisturizer or massage the head.

※Abnormal massage head or posture due to hair or clothing:

Suggestion: the machine should be close to the neck skin and adjust the posture to ensure the correct posture.

※The product is powered by battery, and the voltage is only 3-4V, so please feel free to use it. Suggestion: change massage way, mild from weak to strong, comfortable intensity is appropriate. (Note: pulse massage requires creeping conductive muscle acupoint deep massage to stimulate blood, accelerate circulation, and prick. It stimulates the neck muscles to recover their elasticity. Slight electric flu is not electricity leakage.)

3. After massage, the neck appears red mark

※The pulse wave can stimulate the neck muscle group and the blood circulation of meridians, because of the high activity of epidermal tissue in some people

Local congestion will occur. Suggestion: hot compress towel for 1 minute or rub hot hand for 3 minutes.

※If the effect is not big, wait for 30 minutes, hyperemia will disappear.

4. How to charge (applicable to hx-1680 Series)

※Please use the USB charging cable provided by the instrument to connect the charging power (output power) of mobile phone, power bank or computer. The voltage is 5V, and the current is no more than 1.5A. When charging, the indicator light flashes, and it takes about 3-5 hours to fill up. After full charging, the indicator light is on for a long time

※Do not use this product when charging:

※The service life of the battery will be shortened if the instrument is not used for a long time. It is recommended that the battery be fully charged at least once a month.

5. Why there is no charging adapter for the product (hx-1680 (charging type) series is applicable)

※The product supports universal Android USB cable for charging.

※It supports the market standard universal mobile phone adapter, computer external USB, power bank, vehicle USB charger and power socket.