Mini Fascia Gun




Operating Instruction

An overview of the product



Handle

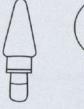



Spherical head





U first

Spike head

Suitable for impacting Applicable to massage Applicable to massage It is suitable for deep tissue, plantar all muscle parts, such the neck and cervical acupoints, such as as the palm of the foot, Achilles tendon. meridian joints, etc


Conical head

massaging large muscle groups, legs, buttocks and other deep stimulation.

Products use

Open/joint (button type)


- Long press the start button for three seconds
- Long press the start button for three seconds to shut down. Efforts to regulate
- The starting gear is in first gear.
- Click the start button once to add
- Maximum gear is in six gear Click the start key once to change to first gear.

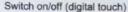

- Long press the "+" button for three seconds to boot
- Long press the "+" button for three seconds to shut down Efforts to regulate
- Power on at the initial gear level 0
- Click "+-" to adjust the strength
- Max. 32 gear




Home key

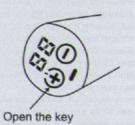


Note: in order to protect the bones and muscles, avoid long-term fixed point massage, resulting in Intelligent protection mode that automatically stops after 15 minutes of damage and closing.

Replace massage head

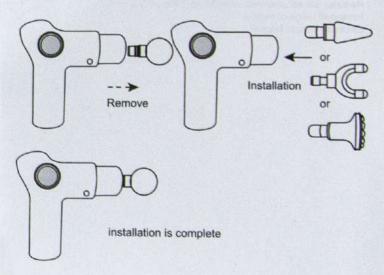

About the charge

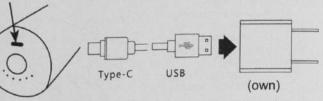

Note: Input is 5V-2A mobile phone charging head

- 1. Power on/off button (long press to power on/long press to power off)
- 2. Press the open key to switch gears
- massage head can be disassembled, replaceable, remove the massage head, as long as the force to pull it out.
- 4. Please use a 5V charging device to charge
- 5. Please charge for three hours for the first time. Please connect the 5V power supply line with the data cable to charge

Button model: 6 lights sparkle when charging, and 6 lights are always on when charging

Touch model: when the product is charged, the screen displays the quantity of electricity and the number flashes. When the product is fully charged, the screen displays "99", which means it is fully charged

Indications

- Pain and spasms from muscle contusions, sprains, strains
- Helps the edema fluid flow to the swollen area
- Loosen and thickened connective tissue and fascia
- Reduces the accumulation of lactic acid in muscles
- Increased range of motion
- Eliminate muscle fatique

Taboo

- Aneurysms, bleeding, using blood thinners
- Heart disease, pacemaker or defibrillator
- Pregnancy, cancer
- The device should not be used within 3 inches of internal fixation for up to 90 days
- Sensitive areas: head, face, cervical vertebra, spinal vertebra
- The tibia, dorsal foot, dorsal hand, etc., where the muscle cover is less near the bone

Listing

Host x1

Type-Cx1

Instructions x1







Spherical head







Spike head

Conical head



1

Product name Fascial gun

certificate

Inspection job No



Product name: Mini fascia gun Rated voltage: 7.4V Motor power: 30W

Interface Type: type-c RPM: 1800-3200 RPM

Gear position: 6 levels adjustment (button type Gear position: 32-gear adjustment (touch type

Massage head number: 4

Material: electronic components +ABS plas

Net weight: 0.46 KG Product size: 142x56x140mm

Warranty Card

Serial Number

Purchase Date

Customer Name

Customer Tel

Customer Email

Warranty Service Terms:

- (1) The warranty of this product is one year, and the built-in battery is guaranteed for six months. (2) The warranty starts from the first purchase date.
- (2) If there is a quality problem that requires warranty, please send this warranty card. Purchase documents and products

Send it back together to the dealer or engineering repair service.

The following conditions are not covered by the warranty:

- (1) The warranty card that cannot pass the formal purchase certificate and fill in the detailed information.
- (2) damage caused by factors considered or irresistible.
- (3) It is believed that the disassembly of the fascia gun caused irreversible Shanghai.

The Company reserves the right to the final interpretation of the above-mentioned after-sales terms and restrictions

