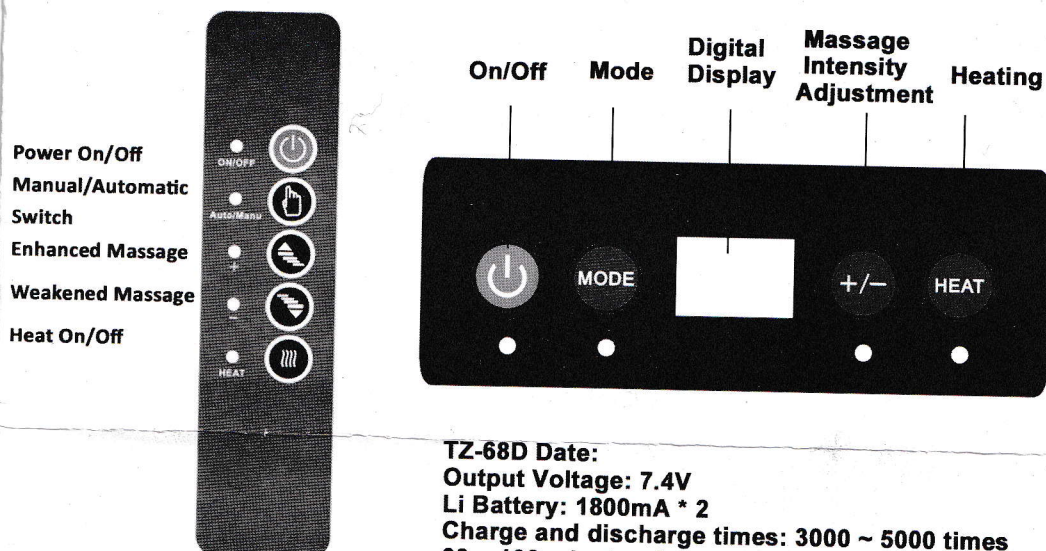




TZ-68D Date:

Output Voltage: 7.4V

Li Battery: 1800mA * 2

**Charge and discharge times: 3000 ~ 5000 times
60 ~ 100 minutes Continuous working time after
full charge**

Net Weight: 0.5kg

Instructions for use:

- 1. Wear:** If the patient's hand with flexion, first open the patient's finger, put the ball into the hands of patients, let the patient naturally bent fingers holding the ball, and fixed fingers.
- 2. Operation of the machine:** Connect to the power, press Power on button,
- 3. Heating Function:** Press the Heat On/Off button, light the indicating light. The infrared light in the ball is turned on. The ball started getting hot. Because the infrared lamp is inside of the ball, the ball's silicone outer thermal slow, usually the ball become warm after 2 minutes.
- 4. Manual and automatic mode switching:** The default mode is automatically. Ball starts vibrating, automatically change from high to low intensity. Press Manual/Automatic Switch, change the machine to Manual mode, Adjustable vibration intensity divided into 1-5 files.

Attention:

- 1. Cleaning,** Strap can be washed. The ball and the remote controller can be used wet wipe (not water)
- 2.** If the long time does not use, please unplug the power
- 3.** Read the manual carefully before using. If abnormal, please stop using the machine immediately.



TZ-68A



TZ-68B



TZ-68C



TZ-68D

Hot compress massage ball:
 The original finger passive rehabilitation training
 Let the movement-disorder fingers move
 Finger passive rehabilitation training
 Infrared heat compress massage
 Adjustable 5 kinds of massage strength

Model: TZ-68A.B.C.D

Scope of application:

Stroke, cerebral hemorrhage, cerebral thrombosis, cerebral palsy, cerebral palsy, traumatic brain injury, cervical vertebra, spinal injury caused by other reasons for hand function abnormal, finger dystonia, hand, finger flexion spasms, hand muscle atrophy, etc.

Long time work cause hand fatigue.