

## Use of guidelines


Gravity handle / meridian scraper  
Enjoy happy and perfect life

### Certificate


Examination clerk: .....	
Date of inspection: .....	
Date of production: .....	
Product Code: .....	

### Warranty card

<b>Warranty card</b> Please keep the card properly.	
<small>Note: the products sold by the company are available from the date of purchase. One year warranty, if the warranty period is over, if you still need to provide maintenance service. We will collect fees according to the actual circumstances.</small>	
<small>The following circumstances do not belong to the scope of free maintenance:</small>	
<small>Failed to produce warranty card or purchase original voucher.</small>	
<small>Warranty cards or Invoices have been altered.</small>	
<small>Product damage caused by unexpected factors.</small>	
<small>Failure to operate according to instructions.</small>	
<small>Product damage caused by dismantling and maintenance without authorization.</small>	
product (Model: _____ number: _____)	customer (Name: _____ Telephone: _____)
Exclusive shop (Name: _____ Telephone: _____)	After sale service center (Name: _____ Telephone: _____)
Salesperson (Signature: _____)	address: _____

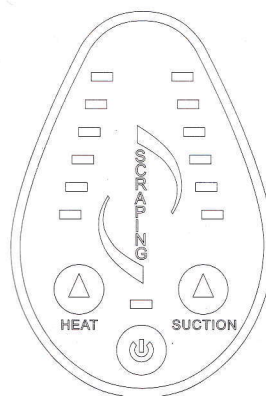

Computer edition

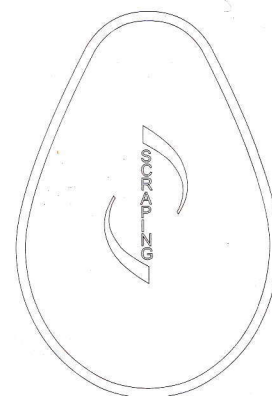

Simple version  
(no keystrokes)

power	Power switch
heating	Temperature 6 adjustable
suction	Suction 6 adjustable
Simplified version can be directly linked to USB without keys.	

## Usage method:

1. clean the skin and apply essential oil.(simplified version skipping the following steps can be directly inserted into USB).
2. plug in the 12V power adapter (equipped with the instrument).
3. when the household power is connected, the power button light is red, and the blue light is pressed.  
At this time, we can adjust the "heat energy" and "negative" according to the needs of customers.Pressure attraction".
4. press the hot energy button and the temperature is divided into 6 files.  
The higher the temperature is, the higher the temperature is.
5. press the button of negative pressure, and the negative pressure is divided into 6 files, which are raised every time.  
The higher the negative pressure is.
6. after using the instrument, press the intermediate power button, then unplug the power supply and clean it.  
After good instrument, place it in dry place.

Note: the default working time of the instrument is 30 minutes, and it will be automatic after 30 minutes.

Suspend work. The instrument will stop working for 20 minutes and then continue to use.

It can extend the life of instrument effectively.

Massage technique

Neck: brush from top to bottom, often scraping and wiping the front part.

The role of disease

Back: scrape and clear from top to bottom.

Lower limbs: dredge up and down

Sole: massage from heel to toe.

Abdomen: the abdomen is dredged up and down, and it can also be operated by rotating circles.

Upper limbs: dredge up and down

Leg: from the shank up scraping wipe dredge, used for leg lymphatic detoxification, shaping leg type.

Eliminate edema. Eliminate leg drying and remove leg fat.

## parameter

Product name	Scraping instrument
Rated voltage	5V or 12V
maximum power	24W
Net weight	≈0.9KG

## Function

1. Negative pressure attraction exercises aerobic exercise on skin and adipose tissue.
2. Scraping function
- 3.continuous heat, relaxation of nerves and release of pressure.
- 4.eliminate fatigue and promote physical recovery.

## Clean maintenance

1. before cleaning the massager, remove the power plug first.
2. do not use benzene, thinner or other corrosive solvents to clean your massager.
3. do not use the wet cloth to clear the power switch of the massager. You can not use alcohol. You must use dry cotton cloth.
4. If you do not use it for a long time, it is recommended that you roll up the power adapter wire and store the massager in a clean and dry environment.
5. do not store the massager in high temperature, open fire or direct sunlight.
6. when using the massager, you should gently handle it gently. Avoid throwing, pressing, bumping or strong impact.
7. It is recommended to replace the filter core after using the scraping instrument for 2-3 times to prevent the functional decline caused by blockage